BODY BEAT


BEDIENUNGSANLEITUNG
USER MANUAL
MODE D'EMPLOI
GUÍA DE USUARIO

SPECIAL EDITION JOEY KELLY

Joey Kelly

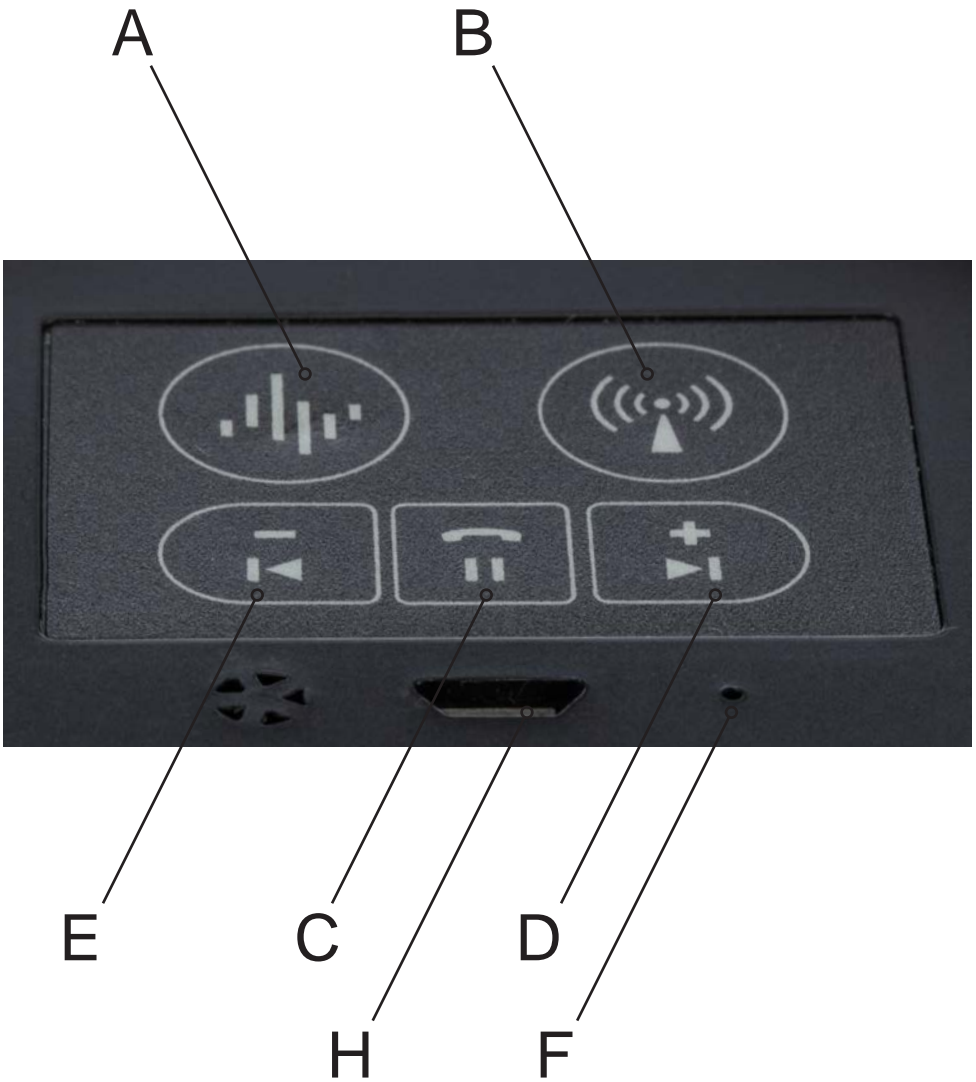WARTUNG & GARANTIE

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MAINTENANCE & WARRANTY

1. Body Beat is not waterproof! Keep away from water and humidity.
2. Clean the Body Beat with a wet cloth only – do not immerse in water and do not use detergents.
3. Avoid dropping or crashing the device as this is not covered by the warranty.
4. Do not expose to direct sunlight!
5. Do not open or disassemble as this invalidates the warranty. In case of quality issues, contact your dealer or the service center.

SPECIFICATIONS

- High intensity vibration, 2 speed setting levels: “Soft” and “Strong”.
- 10 minutes auto-stop function
- Bluetooth® and built-in speakers for playing your favorite music
- Compatible with Micro-USB standard for charging the battery
- If fully charged, device can be used for approx. 3.5 hours (2.5 hours continuous use)
- LED light indicator for set level (“Soft” or “Strong”)
- EVA-foam roller for maximum vibration efficiency
- Acoustic indicator for Bluetooth® function

BODY BEAT BENEFITS

Body Beat is a cutting-edge fitness/rehabilitation device with built-in speakers for daily use.

- Increases the flexibility and motion range by loosening and lengthening the body muscles
- Improves blood-circulation – ideal for warm-up or cool-down before/after physical activity
- Reduces muscle soreness and stiffness and therefore supports the recovery phase
- Effective device for myofascial release

Body Beat can be used as dumbbell for your fitness workout and as vibration foam roller for improving your body`s overall performance. It provides Bluetooth® and speakers for enjoying music and it is suitable for everyday use.

Product dimensions: approx. 308 x 105 x 112 mm

Product weight: approx. 1 kg

Power consumption: Music via Bluetooth®: 0,5 W
Vibration: 3,7 W

User Manual

- A. Press button A once for strong vibration (LED lights up in bright red color). Press button A twice for soft vibration (LED lights up in light red color). A third press on button A turns off the power (LED off). The device will automatically turn off after 10 minutes of usage.
- B. Press button B to start Bluetooth® connection. An acoustic signal is emitted and the LED flashes in red. The device name is “X-BellPro” – after the connection has been established successfully, the LED constantly lights up in red. Press and hold button B for approx. 3 seconds to turn off Bluetooth® connection – an acoustic signal is emitted and the LED will turn off.
- C. Pause button to pause music and answering a phone call.
- D. Press button D to jump to next song or press and hold the button to increase volume.
- E. Press button E to jump to previous song or press and hold the button to decrease volume.
- F. Reset button
- G. Micro-USB port (for charging the battery)
- H. Microphone connection (for phone conversation)

Assembling Instructions

Step 1

Place the foam roll vertically on an even ground (“V” symbol facing up).

Step 2

Slide the dumbbell into the foam roll.

Step 3

Align a corner of the dumbbell with the “V” symbol.

Step 4

Carefully press the dumbbell completely into foam roll using a finger to help the foam perfectly fitting around the end of the dumbbell.

Step 5

The foam roller is completely assembled.

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MAX Trader GmbH